



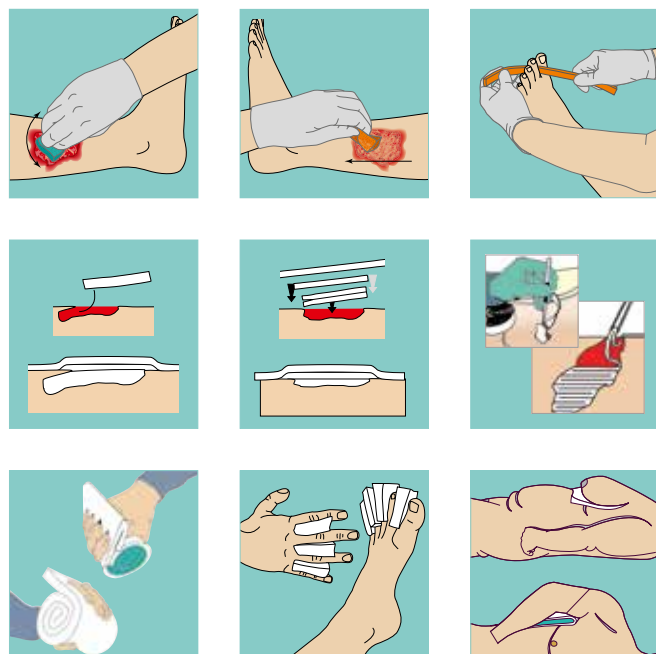
Product List

Edition April 2023

for Pharmacies and Medical Specialised Trade

LIGASANO®

Simply works



Wound Cleaning - Wound Healing - Relapse Prevention

Table of Content:

Wound cleaning	page 3
Sterile wound treatment	page 4 - 5
Wound filler for NPWT	page 6
Unsterile wound treatment	page 7
Large sheets	page 8
Special dressings (e.g. compresses with slit, toe dressing, palm dressing)	page 9 - 11
Further products and special items	page 12 - 14
Appendix	
Brief overview of the products, sorted by article numbers	
Brief overview of the products, sorted by UDI-DI numbers	
Excerpts from the instructions for use	
Terms and conditions	

Please pay attention:

Always read the application instructions enclosed with the original packaging before use! The product descriptions given in this catalog do not release the user from reading the application manual.

The LIGASANO Product Family:

○ = LIGASANO® white
 ● = LIGASANO® orange
 ● = LIGASANO® green

**LIGASANO® white:**

- mixed-pored, comparatively macro-pored polyurethane foam
- cell type 750

Properties:

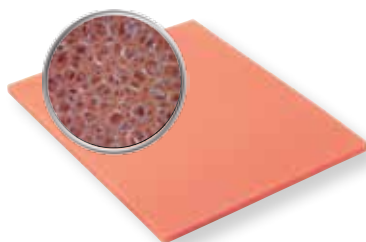
- rapidly falling pressure tension



- mechanical stimulus



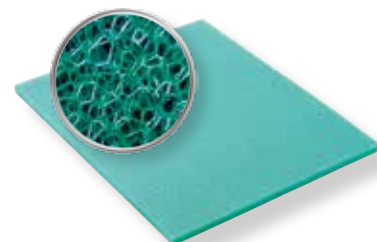
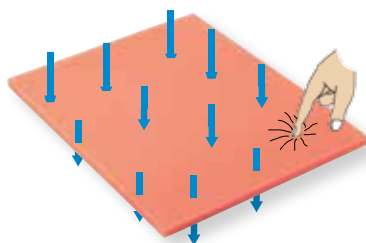
- controlled suction effect

**LIGASANO® orange:**

- open-pored, macropored cell grid with removed cell membranes
- cell type 850

Properties:

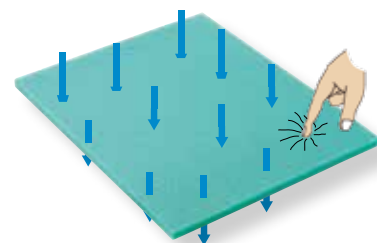
- open-celled and macropored structure
- medium-strong mechanical stimulus
- nearly completely permeable to steam

**LIGASANO® green:**

- open-pored, very macropored cell grid with removed cell membranes
- cell type 1500

Properties:

- open-celled and macropored structure
- intense mechanical stimulus
- nearly completely permeable to steam

**LIGASANO® Sterile Available in Dispenser Boxes:**

Several sterile LIGASANO® articles are available in a functional dispenser box.

The dispensers are perfectly suited for doctor's examination room and the dressing trolley. But also in other situations it facilitates the hygienic handling of LIGASANO®.

The available dispenser boxes are marked in this price list. Please see the boxes in the section of sterile LIGASANO® products for wound cleaning and wound-treatment on page 3-6 and in the section of special dressings on page 10.



Wound Cleaning

LIGASANO® Wundputzer® (Wound Cleaning Sponges) for Mechanical Debridement

Debridement is a key component of wound treatment. The sterile LIGASANO® wound dressings with the additional name „Wundputzer®“ are intended for only a short-term use on or in the wound, to get a clean wound rapidly. A clean wound is one of the preconditions for wound healing.

Depending on kind and amount of the coatings we have three different wound cleaning pads, a combination may make sense. If you use **LIGASANO® white as wound dressing, it causes a mechanical wound cleaning. But persistent coatings are removed more effective and quickly with LIGASANO® green and/or orange.** LIGASANO® Wundputzer® wound cleaning sponges lead to an optimal cleaning result, both in **rotation technique** (180° rotation) and in **wiping technique** (wipe in one direction only).



Wundputzer® intensiv

Removes very stubborn, tough, solid and/or necrotic coatings, can absorb a lot of detritus.



Wundputzer® medium

Removes relatively tough, solid and/or necrotic coatings, can absorb a lot of detritus.



Wundputzer® soft

Removes soft coatings and biofilm.

Item no.	UDI-DI	Shipping unit		
16814-010	111725184144	10 pcs. 5 x 5 x 2 cm green	sterile	●
16814-015	111725185834	15 pcs. 5 x 5 x 2 cm green dispenser box	sterile	●
16814-030	111725186400	30 pcs. 5 x 5 x 2 cm green	sterile	●
16815-010	111725187063	10 pcs. 15 x 5 x 0.6 cm green	sterile	●
36814-010	111725207110	10 pcs. 5 x 5 x 2 cm orange	sterile	●
36814-015	111725208800	15 pcs. 5 x 5 x 2 cm orange dispenser box	sterile	●
36814-030	111725209463	30 pcs. 5 x 5 x 2 cm orange	sterile	●
15342-010	111725126110	10 pcs. 5 x 5 x 2 cm white	sterile	○
15342-015	111725127800	15 pcs. 5 x 5 x 2 cm white dispenser box	sterile	○
15342-030	111725128463	30 pcs. 5 x 5 x 2 cm white	sterile	○
36815-005	111725212513	OP consisting of 5 peel pack each with 1 piece LIGASANO® green 5 x 5 x 2 cm 1 piece LIGASANO® orange 5 x 5 x 2 cm 1 piece LIGASANO® white 5 x 5 x 2 cm	sterile sterile sterile	● ● ○



16812-010	111725182915	10 pcs. 30 x 1.5 x 0.6 cm green Wundputzer® interdigital intensive	sterile	●
16812-028	111725183578	28 pcs. 30 x 1.5 x 0.6 cm green Wundputzer® interdigital intensive dispenser box	sterile	●
36812-010	111725204288	10 pcs. 30 x 1.5 x 0.6 cm orange Wundputzer® interdigital medium	sterile	●
36812-028	111725205978	28 pcs. 30 x 1.5 x 0.6 cm orange Wundputzer® interdigital medium dispenser box	sterile	●
15372-010	111725154561	10 pcs. 50 x 1.5 x 0.4 cm white Wundputzer® interdigital soft	sterile	○
15372-028	111725155127	28 pcs. 50 x 1.5 x 0.4 cm white Wundputzer® interdigital soft dispenser box	sterile	○

LIGASANO® white: Sterile Wound Dressings, Wound Strips, Cavity Dressings

The **rapidly falling pressure tension** of LIGASANO® white reduces the pressure on the wound, promotes granulation and reduces the counterpressure on the new developed granulation tissue. The **mechanical stimulus** (this means that the structure of the dressing is perceived by the tissue) promotes the local blood flow and activates the self-cleaning. Because only surplus exudate is absorbed (**controlled sorption effect, exudate management**), LIGASANO® white does not dry out the wound. Depending of the thickness of the dressing, an intermittent negative pressure evolved. LIGASANO® white can be used for all wound types and wound stages.

Sterile Wound Dressings

Item no.	UDI-DI	Shipping unit
15341-010	111725124978	10 pcs. 5 x 5 x 1 cm sterile ○
15341-040	111725125544	40 pcs. 5 x 5 x 1 cm sterile dispenser box ○
15342-010	111725126110	10 pcs. 5 x 5 x 2 cm sterile ○
15342-015	111725127800	15 pcs. 5 x 5 x 2 cm sterile dispenser box ○
15342-030	111725128463	30 pcs. 5 x 5 x 2 cm sterile ○
15346-010	111725129029	10 pcs. 10 x 10 x 1 cm sterile ○
15346-020	111725130947	20 pcs. 10 x 10 x 1 cm sterile dispenser box ○
15346-200	111725131513	200 pcs. 10 x 10 x 1 cm sterile ○
15347-010	111725132176	10 pcs. 10 x 10 x 2 cm sterile dispenser box ○
15350-015	111725133866	15 pcs. 15 x 10 x 0.5 cm sterile ○
15350-030	111725134432	30 pcs. 15 x 10 x 0.5 cm sterile ○
15351-010	111725135095	10 pcs. 15 x 10 x 1 cm sterile ○
15351-020	111725136785	20 pcs. 15 x 10 x 1 cm sterile ○
15351-140	111725137351	140 pcs. 15 x 10 x 1 cm sterile ○
15356-010	111725145673	10 pcs. 15 x 10 x 2 cm sterile ○
15352-010	111725140401	10 pcs. 24 x 16 x 1 cm sterile ○
15352-070	111725141064	70 pcs. 24 x 16 x 1 cm sterile ○
15357-005	111725146239	5 pcs. 24 x 16 x 2 cm sterile ○
15360-010	111725147929	10 pcs. 29.5 x 24.5 x 1 cm sterile ○
15385-010	111725162883	10 pcs. Ø 5 cm x 1 cm round compress sterile ○
15385-040	111725163449	40 pcs. Ø 5 cm x 1 cm round compress, sterile dispenser box ○



Sterile Bandages and Mini Bandages

Item no.	UDI-DI	Shipping unit
15378-030	111725158046	30 pcs. 30 x 5 x 0.3 cm mini bandage sterile ○
15379-030	111725159736	30 pcs. 50 x 5 x 0.3 cm mini bandage sterile ○
15381-008	111725160530	8 pcs. 300 x 5 x 0.3 cm bandage, width 5 cm sterile ○
15383-004	111725161193	4 pcs. 300 x 10 x 0.3 cm bandage, width 10 cm sterile ○



Sterile Wound Strips / Cavity Dressings

Item no.	UDI-DI	Shipping unit
15370-010	111725148592	10 pcs. 300 x 2.5 x 0.4 cm wound strip sterile ○
15371-007	111725151642	7 pcs. 100 x 1.5 x 0.4 cm wound strip mini sterile ○
15371-014	111725152208	14 pcs. 100 x 1.5 x 0.4 cm wound strip mini, dispenser box sterile ○
15371-028	111725153995	28 pcs. 100 x 1.5 x 0.4 cm wound strip mini sterile ○
15372-010	111725154561	10 pcs. 50 x 1.5 x 0.4 cm wound strip micro sterile ○
15372-028	111725155127	28 pcs. 50 x 1.5 x 0.4 cm wound strip micro, dispenser box sterile ○
15375-010	111725156817	10 pcs. 6 x 2.5 x 0.4 cm sticks sterile ○
15375-060	111725157480	60 pcs. 6 x 2.5 x 0.4 cm sticks sterile ○



LIGASANO® green and orange: Sterile Wound Dressing and Wound Filler for NPWT

The **mechanical stimulus** (this means that the structure of the dressing is perceived by the tissue) promotes the local blood flow and activates the self-cleaning. Because of its open-pored structure LIGASANO® green and LIGASANO® orange **do not show any sorption effect**, i.e. you need a suction source for an active wound drainage or NPWT. By its open-pored and great-pored structure, viscous exudate and detritus is drained very well. LIGASANO® green and LIGASANO® orange are both **almost completely permeable to steam**. In the case of very poor exuding wounds LIGASANO® may stick to the wound by drying exudate.

LIGASANO® green: Sterile Wound Dressings and Wound Filler for NPWT

Item no.	UDI-DI	Shipping unit	
16814-010	111725184144	10 pcs. 5 x 5 x 2 cm	sterile ●
16814-015	111725185834	15 pcs. 5 x 5 x 2 cm dispenser box	sterile ●
16814-030	111725186400	30 pcs. 5 x 5 x 2 cm	sterile ●
16815-010	111725187063	10 pcs. 15 x 10 x 0,6 cm	sterile ●
16816-010	111725188753	10 pcs. 15 x 10 x 1 cm	sterile ●
16816-020	111725189319	20 pcs. 15 x 10 x 1 cm	sterile ●
16817-010	111725191803	10 pcs. 15 x 10 x 2 cm	sterile ●
16820-005	111725193032	5 pcs. 24 x 16 x 2 cm	sterile ●
16820-035	111725194722	35 pcs. 24 x 16 x 2 cm	sterile ●
16812-010	111725182915	10 pcs. 30 x 1,5 x 0,6 cm wound strip mini	sterile ●
16812-028	111725183578	28 pcs. 30 x 1,5 x 0,6 cm wound strip mini, dispenser box	sterile ●
16865-010	111725199994	10 pcs. 200 x 2 x 0,6 cm wound strip	sterile ●



LIGASANO® orange: Sterile Wound Dressings and Wound Filler for NPWT

Item no.	UDI-DI	Shipping unit	
36814-010	111725207110	10 pcs. 5 x 5 x 2 cm	sterile ●
36814-015	111725208800	15 pcs. 5 x 5 x 2 cm dispenser box	sterile ●
36814-030	111725209463	30 pcs. 5 x 5 x 2 cm	sterile ●
36817-010	111725213176	10 pcs. 15 x 10 x 2 cm	sterile ●
36820-005	111725214866	5 pcs. 24 x 16 x 2 cm	sterile ●
36812-010	111725204288	10 pcs. 30 x 1,5 x 0,6 cm wound strip mini	sterile ●
36812-028	111725205978	28 pcs. 30 x 1,5 x 0,6 cm wound strip mini, dispenser box	sterile ●



Unsterile Wound Treatment

LIGASANO® white: Unsterile Dressing for Concomitant Wound Treatment, Secondary Dressing, Cushioning, for Prevention and Therapy of Intertrigo

The **rapidly falling pressure tension** of LIGASANO® white enables almost even distribution of pressure from the surface padding and thus avoids the dreaded pressure points. All LIGASANO® products can simply be cut to fit. This is how aids such as pulleys, wedges and heel shoes, for example, are created. The **mechanical stimulus** of LIGASANO® white promotes the peripheral blood supply where there is contact to skin, and the patients residual pressure increases. Surplus sweat is absorbed, leaving normal skin moisture, and care of the skin is also associated with the mechanical stimulus. The **controlled sorption effect** prevents maceration of the skin effectively. The preservation of the physiological skin milieu is supported.

In bedridden patients and when joints are immobilised, contracture (joint stiffness) and muscle atrophy (shrinkage of the muscles) may develop. Patients with inflammatory or degenerative joint diseases are particularly at risk, but also patients whose joints are immobilised, e.g. due to protective posture, paralysis or weakness. Measures to maintain an intact musculo-skeletal system and functional joints include movement exercises, mobilisation and suitable positioning.

Unsterile Dressing, Smaller Sizes

Item no.	UDI-DI	Shipping unit
15151-010	111725106078	10 pcs. 15 x 15 x 1 cm unsterile
15152-010	111725107768	10 pcs. 24 x 16 x 1 cm unsterile
15153-005	111725108334	5 pcs. 24 x 16 x 2 cm unsterile
15158-013	111725110818	13 pcs. 15 x 10 x 2 cm unsterile
15008-104	111725083015	104 pcs. 15 x 10 x 2 cm unsterile sub-packed 8 x 13 pcs.
15159-026	111725111481	26 pcs. 15 x 10 x 1 cm unsterile
15009-208	111725084705	208 pcs. 15 x 10 x 1 cm unsterile sub-packed 8 x 26 pcs.



Unsterile Bandages

Item no.	UDI-DI	Shipping unit
15571-012	111725171771	12 pcs. 300 x 5 x 0.3 cm unsterile
15573-006	111725174690	6 pcs. 300 x 10 x 0.3 cm unsterile



Unsterile Dressing for Special Applications

Item no.	UDI-DI	Shipping unit
15070-012	111725094256	12 pcs. 6 x 3.5/2.5/1.5 x 0.5 cm unsterile tubes for toes
15075-002	111725095946	2 pcs. Ø 4.5/4.0 x 12 cm unsterile dressing for hallux valgus and big toe
15075-006	111725096512	6 pcs. Ø 4.5/4.0 x 12 cm unsterile dressing for hallux valgus and big toe
15075-024	111725097175	24 pcs. Ø 4.5/4.0 x 12 cm unsterile dressing for hallux valgus and big toe
15078-012	111725098865	12 pcs. Ø 2.5 x 12 cm unsterile dressing for hammer toes
19100-002	111725200706	2 pcs. Ø 15cm x 8/5 cm unsterile ring dressing



LIGASANO® white: Unsterile Dressings for Concomitant Wound Treatment, Large Sheets

Item no.	UDI-DI	Shipping unit	
12001-002	111725062320	2 pcs. 59 x 49 x 2 cm unsterile	○
12001-007	111725064673	7 pcs. 59 x 49 x 2 cm unsterile	○
12001-028	111725065239	28 pcs. 59 x 49 x 2 cm unsterile	○
12002-004	111725066929	4 pcs. 59 x 49 x 1 cm unsterile	○
12002-014	111725067592	14 pcs. 59 x 49 x 1 cm unsterile	○
12002-056	111725068158	56 pcs. 59 x 49 x 1 cm unsterile	○
12003-025	111725069848	25 pcs. 59 x 49 x 0.5 cm unsterile	○
12003-100	111725070642	100 pcs. 59 x 49 x 0.5 cm unsterile	○
12006-001	111725071208	1 pcs. 200 x 100 x 2 cm unsterile	○
12007-002	111725072995	2 pcs. 200 x 100 x 1 cm unsterile	○
13004-005	111725075817	5 pcs. 59 x 98 x 2 cm unsterile	○
13004-020	111725077046	20 pcs. 59 x 98 x 2 cm unsterile	○
13005-010	111725078736	10 pcs. 59 x 98 x 1 cm unsterile	○
13005-040	111725079302	40 pcs. 59 x 98 x 1 cm unsterile	○

**LIGASANO® Starter Unit White + Green**

Item no.	UDI-DI	Shipping unit	
12014-010	111725073561	2 pcs. 59 x 49 x 2 cm unsterile	○
		6 pcs. 59 x 49 x 1 cm unsterile	○
		2 pcs. 55 x 45 x 2 cm unsterile	●
		1 piece zip bag	
12013-003	111767495572	1 piece 59 x 49 x 2 cm unsterile	○
		1 piece 59 x 49 x 1 cm unsterile	○
		1 piece 55 x 45 x 2 cm unsterile	●

**LIGASANO® green: Unsterile Dressings for Concomitant Wound Treatment, Large Sheets**

The surface of LIGASANO® green is relatively rough. If it is used in wounds or on skin lesions may occur, especially by friction. It is advisable to use a thin layer of LIGASANO® white between the skin and LIGASANO® green or an air-permeable textile. Due to its structure, the dressing LIGASANO® green counteracts maceration of the skin effectively. Direct contact to skin is not recommended for this application!

Item no.	UDI-DI	Shipping unit	
16001-002	111725175256	2 pcs. 55 x 45 x 2 cm unsterile	●
16001-007	111725176946	7 pcs. 55 x 45 x 2 cm unsterile	●
16001-028	111725177512	28 pcs. 55 x 45 x 2 cm unsterile	●
16006-001	111725178175	1 pcs. 190 x 90 x 2 cm unsterile	●
16008-006	111725179865	6 pcs. 190 x 45 x 2 cm unsterile operating-table-paddings, single packed	●
16210-002	111725181225	2 pcs. 49 x 49 x 9/1 cm unsterile	●



Palm Dressing

Item no.	UDI-DI	Shipping unit
15090-012	111725101930	10 pcs. 29 x 12 x 1 cm unsterile palm dressing made of LIGASANO® white 2 pcs. 10 cm width, ø 8 cm unsterile spiral rolls made of LIGASANO® green
15091-010	111725102593	10 pcs. 29 x 12 x 1 cm unsterile palm dressing made of LIGASANO® white



Plantar Dressing for Diabetic Feet

Item no.	UDI-DI	Shipping unit
15098-007	111725103159	4 pcs. 35 x 15 x 1 cm unsterile 1 piece 35 x 15 x 2 cm unsterile 1 piece 35 x 15 x 2 cm unsterile 1 piece stocking for fixation



Dressing for Fingers and Toes

Item no.	UDI-DI	Shipping unit
15215-001	111725118912	1 piece 200 x 14 x 0.5 cm unsterile rolled, can cut to fit, with adhesive stripe, for ca. 40-50 tubes for toes



** DBGM Nr. 20 2015 007 500

Sterile Tubes + Caps for Toes

Item no.	UDI-DI	Shipping unit
15391-012	111725164015	12 pcs. 3.5/2.5 x 6 cm sterile
15392-012	111725165705	12 pcs. 2.5/1.5 x 6 cm sterile
15393-012	111725166368	12 pcs. 1.5 x 6 cm sterile
15394-015	111725169287	15 pcs. 14 x 10 x 0.5 cm sterile T-shape dressing for toes, with adhesive stripe for toe diameter from 1.0 to 3.8 cm



** DBGM Nr. 20 2015 007 500

Sterile Dressing for Hallux valgus

Item no.	UDI-DI	Shipping unit
15395-004	111725170081	4 pcs. ø 4.5/4.0 x 12 cm sterile



LIGASANO® white Dressing for Heels, sterile and unsterile

Item no.	UDI-DI	Shipping unit
15330-010	111725123288	10 pcs. Ø 20 x 1 cm sterile ○
15030-010	111725087624	10 pcs. Ø 20 x 1 cm unsterile ○



LIGASANO® Roll Dispenser Boxes without adhesive edge

Item no.	UDI-DI	Shipping unit
15025-001	111725085368	1 pcs. 200 x 2.5 x 1 cm unsterile ○
15225-001	111725120369	1 pcs. 300 x 2.5 x 0.6 cm unsterile ○
15050-001	111725088287	1 pcs. 300 x 5.0 x 0.3 cm unsterile ○
15055-001	111725089977	1 pcs. 200 x 5.5 x 1 cm unsterile ○
15255-001	111725122625	1 pcs. 300 x 5.5 x 0.6 cm unsterile ○
15110-001	111725105415	1 pcs. 200 x 11 x 1 cm unsterile ○
15210-001	111725117222	1 pcs. 300 x 11 x 0.6 cm unsterile ○
16110-001	111725180659	1 pcs. 200 x 11 x 1 cm unsterile ○

with adhesive edge

15215-001	111725118912	1 pcs. 200 x 14 x 0.5 cm unsterile ○
15216-001	111725119575	1 pcs. 200 x 4,5 x 0.4 cm unsterile ○



LIGASANO® white Spiral Dressings

Item no.	UDI-DI	Shipping unit
15217-010	18044765	10 pcs. LIGASANO® white unsterile spiral dressing, Ø 1.6 cm, length 12 cm ○
15218-009	18044759	9 pcs. LIGASANO® white unsterile spiral dressing, Ø 2.5 cm, length 12 cm ○



LIGASANO® Compresses with Slit

Item no.	UDI-DI	Shipping unit
15164-020	111725112047	20 pcs. 5 x 5 x 1 cm sterile ○
15165-020	111725113737	20 pcs. 7.5 x 7.5 x 1 cm sterile dispenser box ○
15166-020	111725116656	20 pcs. 10 x 10 x 1 cm sterile dispenser box ○
15063-020	111725090771	20 pcs. 5 x 5 x 1 cm unsterile ○
15064-020	111725091337	20 pcs. 7.5 x 7.5 x 1 cm unsterile ○
15065-020	111725093690	20 pcs. 10 x 10 x 1 cm unsterile ○



Special Dressings

LIGASANO® white Ear Dressing, unsterile

Item no.	UDI-DI	Shipping unit
15072-018	111752685771	18 pcs. 20 x 10 x 1 cm unsterile ear dressing with punched opening for the ear, for folding <input type="radio"/>



Omega Upholstery Dressing

Item no.	UDI-DI	Shipping unit
36050-006	111725203625	1 piece LIGASANO® orange unsterile Omega dressing, adjustable in 5 heights <input checked="" type="radio"/> 5 pcs. LIGASANO® white unsterile skin contact dressing 15 x 10 x 1 cm <input type="radio"/>

* DBGM Nr. 20 2017 000 091



Trapeze Upholstery Dressing

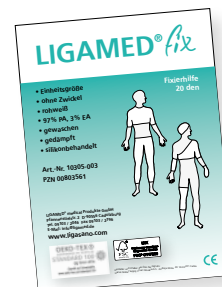
Item no.	UDI-DI	Shipping unit
36060-006	111767494909	1 piece LIGASANO® orange unsterile Trapeze dressing, adjustable in 4 heights <input checked="" type="radio"/> 5 pcs. LIGASANO® white unsterile skin contact dressing 30 x 12 x 1 cm <input type="radio"/>



LIGAMED® *fix* Fixation Aid

Item no.	UDI-DI	Shipping unit
10305-003	110080356115	3 pcs. fixation aid

LIGAMED® *fix* is a practical, glue-free fixation for LIGASANO® for the whole body. It offers high wearing comfort, is easy to adapt and often reusable for several times. Thus you can avoid many applications with adhesive plaster and accordingly its disadvantages.



LIGAMED® Heel Shoe

Item no.	UDI-DI	Shipping unit
32210-001	110467815388	heel shoe, one piece, one size, for all feet, washable, sterilisable; 27,5 x 9 x 18 cm
32211-001	111638996842	heel shoe plus, one piece, one size, for all feet, washable, sterilisable; 27,5 x 9 x 23 cm



Heel Shoe
Item 32210-001

Heel Shoe Plus
Item 32211-001

LIGASANO® Zip Bag

Item no.	UDI-DI	Shipping unit
12000-001	4031463120000	1 piece zip bag 60 x 50 x 14 cm

For hygienic and dust-free storage of LIGASANO® sheets.



LIGAMED® Medical Instruments

Item no.	UDI-DI	Shipping unit
10433-001	4031463104338	1 piece bandage scissors according to Lister, 18 cm, unsterile, certified stainless quality
10441-001	110919242679	1 piece forceps, anatomic, 14.5 cm, unsterile, certified stainless quality



Custom-Made Articles

The following articles are not always in stock, but are produced after the order has been received. If you are interested, please ask for the current delivery time.

Item no.	UDI-DI	Shipping unit
15370-060	111725149158	60 pcs. 300 x 2.5 x 0.4 cm sterile ○
15083-025	111725099431	25 pcs. Ø 3.5 cm x 1 cm unsterile round compress ○
15085-036	111725100240	36 pcs. Ø 5 cm x 1 cm unsterile round compress ○
15099-020	111725104849	20 pcs. 35 x 15 x 1 cm unsterile Plantar dressing for diabetics ○
19100-024	111725201369	24 pcs. Ø 15 cm x 8/5 cm unsterile ring ○
12003-100	111725070642	100 pcs. 59 x 49 x 0.5 cm unsterile ○
13008-006	111725082449	6 pcs. 200 x 50 x 2 cm unsterile OP table sheets, single packed ○
16817-070	111725192466	70 pcs. 15 x 10 x 2 cm sterile ●
16008-006	111725179865	6 pcs. 190 x 45 x 2 cm unsterile OP table sheets, single packed ●
36120-004	4031463361205	4 pcs. 12 x 12 x 12 cm unsterile ●
36121-002	4031463361212	2 pcs. 24 x 12 x 12 cm unsterile ●



LIGASANO® Ambulance Package

Item no.	UDI-DI	Shipping unit
15353-015	110595893796	1 piece 24 x 16 x 1 cm sterile ○
		1 piece 15 x 10 x 1 cm sterile ○
		1 piece 15 x 10 x 0.5 cm sterile ○
		1 piece 10 x 10 x 1 cm sterile ○
		1 piece 5 x 5 x 1 cm sterile ○
		1 piece 6 x 2.5 x 0.4 cm sterile ○
		1 piece 300 x 2.5 x 0.4 cm sterile ○
		1 piece 100 x 1.5 x 0.4 cm sterile ○
		1x3 pcs. 5 x 5 x 2 cm sterile ●
		1 piece Ø4.5/4.0 x 12 cm sterile ○
		1 piece T-formed 14 x 10 cm sterile ○
		1 piece 300 x 5 x 0.3 cm unsterile ○
		1 piece 300 x 10 x 0.3 cm unsterile ○
		1 piece 200 x 2.5 x 1 cm unsterile ○
		1 piece 6 x 3,5/2,5/1,5 x 1 cm unsterile ○



LIGASANO® Ambulance Package for Diabetologists and Podiatrists

Item no.	UDI-DI	Shipping unit					
15354-016	111725142754	1 piece	5 x	5 x	1 cm	sterile	<input type="radio"/>
		1 piece	10 x	10 x	1 cm	sterile	<input type="radio"/>
		1 piece	15 x	10 x	0.5 cm	sterile	<input type="radio"/>
		1 piece	100 x	1.5 x	0.4 cm	sterile	<input type="radio"/>
		1 piece	6 x	2.5 x	0.4 cm	sterile	<input type="radio"/>
		1 piece	Ø 5.0 x		1 cm	sterile	<input type="radio"/>
		1 piece	Ø 3.5/2.5 x		6 cm	sterile	<input type="radio"/>
		1 piece	Ø 2.5/1.5 x		6 cm	sterile	<input type="radio"/>
		1 piece	Ø 1.5 x		6 cm	sterile	<input type="radio"/>
		1 piece	5 x	5 x	1 cm	unsterile	<input type="radio"/>
		1 piece	10 x	10 x	1 cm	unsterile	<input type="radio"/>
		1 piece	Ø 3.5/2.5/1.5 x		6 cm	unsterile	<input type="radio"/>
		1 piece	Ø 4.5/4.0 x		12 cm	unsterile	<input type="radio"/>
		1 piece	Ø 3.5 x		1 cm	unsterile	<input type="radio"/>
		1 piece	Ø 5.0 x		1 cm	unsterile	<input type="radio"/>
		1 piece	300 x	5 x	0.3 cm	unsterile	<input type="radio"/>



LIGASANO® Universal Package for Podiatrists

12016-007	111725074127	3 pcs.	49 x 29.5 x	1 cm	unsterile	<input type="radio"/>
		2 pcs.	55 x 27.5 x	1 cm	unsterile	<input checked="" type="radio"/>



LIGASANO® Mixed Package for Podiatrists

12017-022	111767496138	5 pcs.	5 x	5 x	1 cm	unsterile	<input type="radio"/>
		5 pcs.	15 x	10 x	1 cm	unsterile	<input type="radio"/>
		5 pcs.	Ø 3.5cm x		1 cm	unsterile	<input type="radio"/>
		1 pcs.	300 x	5 x	0.3 cm	unsterile	<input type="radio"/>
		4 pcs.	Ø 3.5/2.5/1.5 x		6 cm	unsterile	<input type="radio"/>
		2 pcs.	Ø 4.5/4.0 x		12 cm	unsterile	<input type="radio"/>



Short Overview (Numeric Sorting)

Art. no.	UDI-DI	Description	Size	Pcs.
10305-005	110080356115	LIGAMED fix	XL	3
10433-001		Lister scissors	18 cm	1
10441-001	110919242679	Forceps, anatomic	14.5 cm	1
12000-001		Zip Bag	60 x 50 x 14 cm	1
12001-002	111725062320	LIGASANO® white unsterile	59 x 49 x 2 cm	2
12001-007	111725064673	LIGASANO® white unsterile	59 x 49 x 2 cm	7
12001-028	111725065239	LIGASANO® white unsterile	59 x 49 x 2 cm	28
12002-004	111725066929	LIGASANO® white unsterile	59 x 49 x 1 cm	4
12002-014	111725067592	LIGASANO® white unsterile	59 x 49 x 1 cm	14
12002-056	111725068158	LIGASANO® white unsterile	59 x 49 x 1 cm	56
12003-025	111725069848	LIGASANO® white unsterile	59 x 49 x 0.5 cm	25
12003-100	111725070642	LIGASANO® white unsterile	59 x 49 x 0.5 cm	100
12006-001	111725071208	LIGASANO® white unsterile	200 x 100 x 2 cm	1
12007-002	111725072995	LIGASANO® white unsterile	200 x 100 x 1 cm	2
12013-003	111767495572	LIGASANO® Combi pack	Various sizes	3
12014-010	111725073561	LIGASANO® Combi pack	Various sizes	10
12016-007	111725074127	LIGASANO® Combi pack	Various sizes	7
12017-022	111767496138	LIGASANO® Combi pack	Various sizes	22
13004-005	111725075817	LIGASANO® white unsterile	59 x 98 x 2 cm	5
13004-020	111725077046	LIGASANO® white unsterile	59 x 98 x 2 cm	20
13005-010	111725078736	LIGASANO® white unsterile	59 x 98 x 1 cm	10
13005-040	111725079302	LIGASANO® white unsterile	59 x 98 x 1 cm	40
13008-006	111725082449	LIGASANO® white unsterile	200 x 50 x 2 cm	6
15008-104	111725083015	LIGASANO® white unsterile	15 x 10 x 2 cm	104
15009-208	111725084705	LIGASANO® white unsterile	15 x 10 x 1 cm	208
15025-001	111725085368	LIGASANO® white unsterile	200 x 2.5 x 1 cm	1
15030-010	111725087624	LIGASANO® white unsterile	Ø 20 x 1 cm	10
15041-000	110163294359	LIGASANO® white unsterile	5 x 5 x 1 cm	1680
15046-420	110163293793	LIGASANO® white unsterile	10 x 10 x 1 cm	420
15050-001	111725088287	LIGASANO® white unsterile	300 x 5 x 0.3 cm	1
15055-001	111725089977	LIGASANO® white unsterile	200 x 5.5 x 1 cm	1
15063-020	111725090771	LIGASANO® white ust. slit	5 x 5 x 1 cm	20
15064-020	111725091337	LIGASANO® white ust. slit	7.5 x 7.5 x 1 cm	20
15065-020	111725093690	LIGASANO® white ust. slit	10 x 10 x 1 cm	20
15070-012	111725094256	LIGASANO® white unsterile	Ø3.5/2.5/1.5 x 6 cm	12
15072-018	1117252685871	LIGASANO® white unsterile	20 x 10 x 1 cm Ohr	18
15075-002	111725095946	LIGASANO® white unsterile	Ø 4.5/4.0 x 12 cm	2
15075-006	111725096512	LIGASANO® white unsterile	Ø 4.5/4.0 x 12 cm	6
15075-024	111725097175	LIGASANO® white unsterile	Ø 4.5/4.0 x 12 cm	24
15078-012	111725098865	LIGASANO® white unsterile	Ø 2.5 x 12 cm	12
15083-025	111725099431	LIGASANO® white unsterile	Ø 3.5 x 1 cm	25
15085-036	111725100240	LIGASANO® white unsterile	Ø 5.0 x 1 cm	36
15090-012	111725101930	LIGASANO® Combi pack	Diverse Größen	12
15091-010	111725102593	LIGASANO® white unsterile	29 x 12 x 1 cm	10
15098-007	111725103159	LIGASANO® Combi pack	Diverse Größen	7
15099-020	111725104849	LIGASANO® white unsterile	35 x 15 x 1 cm	20
15110-001	111725105415	LIGASANO® white unsterile	200 x 11 x 1 cm	1
15151-010	111725106078	LIGASANO® white unsterile	15 x 15 x 1 cm	10
15152-010	111725107768	LIGASANO® white unsterile	24 x 16 x 1 cm	10
15153-005	111725108334	LIGASANO® white unsterile	24 x 16 x 2 cm	5
15158-013	111725110818	LIGASANO® white unsterile	15 x 10 x 2 cm	13
15159-026	111725111481	LIGASANO® white unsterile	15 x 10 x 1 cm	26
15164-020	111725112047	LIGASANO® white st. slit	5 x 5 x 1 cm	20
15165-020	111725113737	LIGASANO® white st. slit	7.5 x 7.5 x 1 cm	20
15166-020	111725116656	LIGASANO® white st. slit	10 x 10 x 1 cm	20
15210-001	111725117222	LIGASANO® white unsterile	300 x 11 x 0.6 cm	1
15215-001	111725118912	LIGASANO® white unsterile	200 x 14 x 0.5 cm	1
15216-001	111725119575	LIGASANO® white unsterile	200 x 4.5 x 0.4 cm	1
15217-010	111804476512	LIGASANO® white unsterile	Ø 1.6 x 12 cm	10
15218-009	111804475946	LIGASANO® white unsterile	Ø 2.5 x 12 cm	9
15225-001	111725120369	LIGASANO® white unsterile	300 x 2.5 x 0.6 cm	1
15255-001	111725122625	LIGASANO® white unsterile	300 x 5.5 x 0.6 cm	1
15330-010	111725123288	LIGASANO® white sterile	Ø 20 x 1 cm	10
15341-010	111725124978	LIGASANO® white sterile	5 x 5 x 1 cm	10
15341-040	111725125544	LIGASANO® white sterile	5 x 5 x 1 cm	40
15342-010	111725126110	LIGASANO® white sterile	5 x 5 x 2 cm	10
15342-015	111725127800	LIGASANO® white sterile	5 x 5 x 2 cm	15
15342-030	111725128463	LIGASANO® white sterile	5 x 5 x 2 cm	30
15346-010	111725129029	LIGASANO® white sterile	10 x 10 x 1 cm	10
15346-020	111725130947	LIGASANO® white sterile	10 x 10 x 1 cm	20
15346-200	111725131513	LIGASANO® white sterile	10 x 10 x 1 cm	200
15347-010	111725132176	LIGASANO® white sterile	10 x 10 x 2 cm	10
15350-015	111725133866	LIGASANO® white sterile	15 x 10 x 0.5 cm	15
15350-030	111725134432	LIGASANO® white sterile	15 x 10 x 0.5 cm	30

Art. no.	UDI-DI	Description	Size	Pcs.
15351-010	111725135095	LIGASANO® white sterile	15 x 10 x 1 cm	10
15351-020	111725136785	LIGASANO® white sterile	15 x 10 x 1 cm	20
15351-140	111725137351	LIGASANO® white sterile	15 x 10 x 1 cm	140
15352-010	111725140401	LIGASANO® white sterile	24 x 16 x 1 cm	10
15352-070	111725141064	LIGASANO® white sterile	24 x 16 x 1 cm	70
15353-015	110595893796	LIGASANO® Combi pack	Various sizes	15
15354-016	111725142754	LIGASANO® Combi pack	Various sizes	16
15356-010	111725145673	LIGASANO® white sterile	15 x 10 x 2 cm	10
15357-005	111725146239	LIGASANO® white sterile	24 x 16 x 2 cm	5
15360-010	111725147929	LIGASANO® white sterile	29.5 x 24.5 x 1 cm	10
15370-010	111725148592	LIGASANO® white sterile	300 x 2.5 x 0.4 cm	10
15370-060	111725149158	LIGASANO® white sterile	300 x 2.5 x 0.4 cm	60
15371-007	111725151642	LIGASANO® white sterile	100 x 1.5 x 0.4 cm	7
15371-014	111725152208	LIGASANO® white sterile	100 x 1.5 x 0.4 cm	14
15371-028	111725153995	LIGASANO® white sterile	100 x 1.5 x 0.4 cm	28
15372-010	111725154561	LIGASANO® white sterile	50 x 1.5 x 0.4 cm	10
15372-028	111725155127	LIGASANO® white sterile	50 x 1.5 x 0.4 cm	28
15375-010	111725156817	LIGASANO® white sterile	6 x 2.5 x 0.4 cm	10
15375-060	111725157480	LIGASANO® white sterile	6 x 2.5 x 0.4 cm	60
15378-030	111725158046	LIGASANO® white sterile	30 x 5 x 0.3 cm	30
15379-030	111725159736	LIGASANO® white sterile	50 x 5 x 0.3 cm	30
15381-008	111725160530	LIGASANO® white sterile	300 x 5 x 0.3 cm	8
15383-004	111725161193	LIGASANO® white sterile	300 x 10 x 0.3 cm	4
15385-010	111725162883	LIGASANO® white sterile	Ø 5 x 1 cm	10
15385-040	111725163449	LIGASANO® white sterile	Ø 5 x 1 cm	28
15391-012	111725164015	LIGASANO® white sterile	Ø 3.5 x 6 cm	12
15392-012	111725165705	LIGASANO® white sterile	Ø 2.5 x 6 cm	12
15393-012	111725166368	LIGASANO® white sterile	Ø 1.5 x 6 cm	12
15394-015	111725169287	LIGASANO® white sterile	T-Form 14 x 10 cm	15
15395-004	111725170081	LIGASANO® white sterile	Ø 4.5/4.0 x 12 cm	4
15571-012	111725171771	LIGASANO® white unsterile	300 x 5 x 0.3 cm	12
15573-006	111725174690	LIGASANO® white unsterile	300 x 10 x 0.3 cm	6
16001-002	111725175256	LIGASANO® green unsterile	55 x 45 x 2 cm	2
16001-007	111725176946	LIGASANO® green unsterile	55 x 45 x 2 cm	7
16001-028	111725177512	LIGASANO® green unsterile	55 x 45 x 2 cm	28
16006-001	111725178175	LIGASANO® green unsterile	190 x 90 x 2 cm	1
16008-006	111725179865	LIGASANO® green unsterile	190 x 45 x 2 cm	6
16110-001	111725180659	LIGASANO® green unsterile	200 x 11 x 1 cm	1
16210-002	111725181225	LIGASANO® green unsterile	49 x 49 x 9/1 cm	2
16812-010	111725182915	LIGASANO® green sterile	30 x 1.5 x 0.6 cm	10
16812-028	111725183578	LIGASANO® green sterile	30 x 1.5 x 0.6 cm	28
16814-010	111725184144	LIGASANO® green sterile	5 x 5 x 2 cm	10
16814-015	111725185834	LIGASANO® green sterile	5 x 5 x 2 cm	15
16814-030	111725186400	LIGASANO® green sterile	5 x 5 x 2 cm	30
16815-010	111725187063	LIGASANO® green sterile	15 x 5 x 0.6 cm	10
16816-010	111725188753	LIGASANO® green sterile	15 x 10 x 1 cm	10
16816-020	111725189319	LIGASANO® green sterile	15 x 10 x 1 cm	20
16817-010	111725191803	LIGASANO® green sterile	15 x 10 x 2 cm	10
16817-070	111725192466	LIGASANO® green sterile	15 x 10 x 2 cm	70
16820-005	111725193032	LIGASANO® green sterile	24 x 16 x 2 cm	5
16820-035	111725194722	LIGASANO® green sterile	24 x 16 x 2 cm	35
16865-010	111725199994	LIGASANO® green sterile	200 x 2 x 0.6 cm	10
19100-002	111725200706	LIGASANO® white unsterile	Ø 15 x 8/5 cm	2
19100-024	111725201369	LIGASANO® white unsterile	Ø 15 x 8/5 cm	24
32210-001	110467815388	LIGAMED® Heel shoe	One size	1
32211-001	111638996842	LIGAMED® Heel shoe Plus	One size	1
36050-006	111725203625	LIGASANO® Combi pack	5-fold adjustable	6
36060-006	111767494909	LIGASANO® Combi pack	4-fold adjustable	6
36120-004		LIGASANO® orange unsterile	12 x 12 x 12 cm	4
36121-002		LIGASANO® orange unsterile	24 x 12 x 12 cm	2
36812-010	111725204288	LIGASANO® orange sterile	30 x 1.5 x 0.6 cm	10
36812-028	111725205978	LIGASANO® orange sterile	30 x 1.5 x 0.6 cm	28
36814-010	111725207110	LIGASANO® orange sterile	5 x 5 x 2 cm	10
36814-015	111725208800	LIGASANO® orange sterile	5 x 5 x 2 cm	15
36814-030	111725209463	LIGASANO® orange sterile	5 x 5 x 2 cm	30
36815-005	111725212513	LIGASANO® Combi pack	5 x 5 x 2 cm	5
36817-010	111725213176	LIGASANO® orange sterile	15 x 10 x 2 cm	10
36820-005	111725214866	LIGASANO® orange sterile	24 x 16 x 2 cm	5

Short Overview (UDI-DI Sorting)

UDI-DI	Art. no.	Description	Size	Pcs.
110080356115	10305-003	LIGAMED® fix	XL	3
110163293793	15046-420	LIGASANO® white unsterile	10 x 10 x 1 cm	420
110163294359	15041-000	LIGASANO® white unsterile	5 x 5 x 1 cm	1680
110467815388	32210-001	LIGAMED® Heel shoe	One size	1
110595893796	15353-015	LIGASANO® Combi pack	Various sizes	15
110919242679	10441-001	LIGAMED® Forceps	14,5 cm	1
111638996842	32211-001	LIGAMED® Heel shoe plus	One size	1
111725062320	12001-002	LIGASANO® white unsterile	59 x 49 x 2 cm	2
111725064673	12001-007	LIGASANO® white unsterile	59 x 49 x 2 cm	7
111725065239	12001-028	LIGASANO® white unsterile	59 x 49 x 2 cm	28
111725066929	12002-004	LIGASANO® white unsterile	59 x 49 x 1 cm	4
111725067592	12002-014	LIGASANO® white unsterile	59 x 49 x 1 cm	14
111725068158	12002-056	LIGASANO® white unsterile	59 x 49 x 1 cm	56
111725069848	12003-025	LIGASANO® white unsterile	59 x 49 x 0.5 cm	25
111725070642	12003-100	LIGASANO® white unsterile	59 x 49 x 0.5 cm	100
111725071208	12006-001	LIGASANO® white unsterile	200 x 100 x 2 cm	1
111725072995	12007-002	LIGASANO® white unsterile	200 x 100 x 1 cm	2
111725073561	12014-010	LIGASANO® Combi pack	Various sizes	10
111725074127	12016-007	LIGASANO® Combi pack	Various sizes	7
111725075817	13004-005	LIGASANO® white unsterile	59 x 98 x 2 cm	5
111725077046	13004-020	LIGASANO® white unsterile	59 x 98 x 2 cm	20
111725078736	13005-010	LIGASANO® white unsterile	59 x 98 x 1 cm	10
111725079302	13005-040	LIGASANO® white unsterile	59 x 98 x 1 cm	40
111725082449	13008-006	LIGASANO® white unsterile	200 x 50 x 2 cm	6
111725083015	15008-104	LIGASANO® white unsterile	15 x 10 x 2 cm	104
111725084705	15009-208	LIGASANO® white unsterile	15 x 10 x 1 cm	208
111725085368	15025-001	LIGASANO® white unsterile	200 x 2.5 x 1 cm	1
111725087624	15030-010	LIGASANO® white unsterile	Ø 20 x 1 cm	10
111725088287	15050-001	LIGASANO® white unsterile	300 x 5 x 0.3 cm	1
111725089977	15055-001	LIGASANO® white unsterile	200 x 5.5 x 1 cm	1
111725090771	15063-020	LIGASANO® white ust. slot	5 x 5 x 1 cm	20
111725091337	15064-020	LIGASANO® white ust. slot	7.5 x 7.5 x 1 cm	20
111725093690	15065-020	LIGASANO® white ust. slot	10 x 10 x 1 cm	20
111725094256	15070-012	LIGASANO® white unsterile	Ø 3.5/2.5/1.5x6 cm	12
111725095946	15075-002	LIGASANO® white unsterile	Ø 4.5/4.0 x 12 cm	2
111725096512	15075-006	LIGASANO® white unsterile	Ø 4.5/4.0 x 12 cm	6
111725097175	15075-024	LIGASANO® white unsterile	Ø 4.5/4.0 x 12 cm	24
111725098865	15078-012	LIGASANO® white unsterile	Ø 2.5 x 12 cm	12
111725099431	15083-025	LIGASANO® white unsterile	Ø 3.5 x 1 cm	25
111725100240	15085-036	LIGASANO® white unsterile	Ø 5.0 x 1 cm	36
111725101930	15090-012	LIGASANO® Combi pack	Various sizes	12
111725102593	15091-010	LIGASANO® white unsterile	29 x 12 x 1 cm	10
111725103159	15098-007	LIGASANO® Combi pack	Various sizes	7
111725104849	15099-020	LIGASANO® white unsterile	35 x 15 x 1 cm	20
111725105415	15110-001	LIGASANO® white unsterile	200 x 11 x 1 cm	1
111725106078	15151-010	LIGASANO® white unsterile	15 x 15 x 1 cm	10
111725107768	15152-010	LIGASANO® white unsterile	24 x 16 x 1 cm	10
111725108334	15153-005	LIGASANO® white unsterile	24 x 16 x 2 cm	5
111725110818	15158-013	LIGASANO® white unsterile	15 x 10 x 2 cm	13
111725111481	15159-026	LIGASANO® white unsterile	15 x 10 x 1 cm	26
111725112047	15164-020	LIGASANO® white st. slot	5 x 5 x 1 cm	20
111725113737	15165-020	LIGASANO® white st. slot	7.5 x 7.5 x 1 cm	20
111725116656	15166-020	LIGASANO® white st. slot	10 x 10 x 1 cm	20
111725117222	15210-001	LIGASANO® white unsterile	300 x 11 x 0.6 cm	1
111725118912	15215-001	LIGASANO® white unsterile	200 x 14 x 0.5 cm	1
111725119575	15216-001	LIGASANO® white unsterile	200 x 4.5 x 0.4 cm	1
111725120369	15225-001	LIGASANO® white unsterile	300 x 2.5 x 0.6 cm	1
111725122625	15255-001	LIGASANO® white unsterile	300 x 5.5 x 0.6 cm	1
111725123288	15330-010	LIGASANO® white sterile	Ø 20 x 1 cm	10
111725124978	15341-010	LIGASANO® white sterile	5 x 5 x 1 cm	10
111725125544	15341-040	LIGASANO® white sterile	5 x 5 x 1 cm	40
111725126110	15342-010	LIGASANO® white sterile	5 x 5 x 2 cm	10
111725127800	15342-015	LIGASANO® white sterile	5 x 5 x 2 cm	15
111725128463	15342-030	LIGASANO® white sterile	5 x 5 x 2 cm	30
111725129029	15346-010	LIGASANO® white sterile	10 x 10 x 1 cm	10
111725130947	15346-020	LIGASANO® white sterile	10 x 10 x 1 cm	20
111725131513	15346-200	LIGASANO® white sterile	10 x 10 x 1 cm	200
111725132176	15347-010	LIGASANO® white sterile	10 x 10 x 2 cm	10
111725133866	15350-015	LIGASANO® white sterile	15 x 10 x 0.5 cm	15
111725134432	15350-030	LIGASANO® white sterile	15 x 10 x 0.5 cm	30
111725135095	15351-010	LIGASANO® white sterile	15 x 10 x 1 cm	10
111725136785	15351-020	LIGASANO® white sterile	15 x 10 x 1 cm	20
111725137351	15351-140	LIGASANO® white sterile	15 x 10 x 1 cm	140
111725140401	15352-010	LIGASANO® white sterile	24 x 16 x 1 cm	10

UDI-DI	Art. no.	Description	Size	Pcs.
111725141064	15352-070	LIGASANO® white sterile	24 x 16 x 1 cm	70
111725142754	15354-016	LIGASANO® white sterile	Various sizes	16
111725145673	15356-010	LIGASANO® white sterile	15 x 10 x 2 cm	10
111725146239	15357-005	LIGASANO® white sterile	24 x 16 x 2 cm	5
111725147929	15360-010	LIGASANO® white sterile	29.5 x 24.5 x 1 cm	10
111725148592	15370-010	LIGASANO® white sterile	300 x 2.5 x 0.4 cm	10
111725149158	15370-060	LIGASANO® white sterile	300 x 2.5 x 0.4 cm	60
111725151642	15371-007	LIGASANO® white sterile	100 x 1.5 x 0.4 cm	7
111725152208	15371-014	LIGASANO® white sterile	100 x 1.5 x 0.4 cm	14
111725153995	15371-028	LIGASANO® white sterile	100 x 1.5 x 0.4 cm	28
111725154561	15372-010	LIGASANO® white sterile	50 x 1.5 x 0.4 cm	10
111725155127	15372-028	LIGASANO® white sterile	50 x 1.5 x 0.4 cm	28
111725156817	15375-010	LIGASANO® white sterile	6 x 2.5 x 0.4 cm	10
111725157480	15375-060	LIGASANO® white sterile	6 x 2.5 x 0.4 cm	60
111725158046	15378-030	LIGASANO® white sterile	30 x 5 x 0.3 cm	30
111725159736	15379-030	LIGASANO® white sterile	50 x 5 x 0.3 cm	30
111725160530	15381-008	LIGASANO® white sterile	300 x 5 x 0.3 cm	8
111725161193	15383-004	LIGASANO® white sterile	300 x 10 x 0.3 cm	4
111725162883	15385-010	LIGASANO® white sterile	Ø 5.0 x 1 cm	10
111725163449	15385-040	LIGASANO® white sterile	Ø 5.0 x 1 cm	40
111725164015	15391-012	LIGASANO® white sterile	Ø 3.5/2.5 x 1 cm	12
111725165705	15392-012	LIGASANO® white sterile	Ø 2.5/1.5 x 1 cm	12
111725166368	15393-012	LIGASANO® white sterile	Ø 1.5 x 1 cm	12
111725169287	15394-015	LIGASANO® white sterile	T, 14 x 10 x 0.5 cm	15
111725170081	15395-004	LIGASANO® white sterile	Ø 4.5/4.0 x 12 cm	10
111725171771	15571-012	LIGASANO® white unsterile	300 x 5 x 0.3 cm	12
111725174690	15573-006	LIGASANO® white unsterile	300 x 10 x 0.3 cm	6
111725175256	16001-002	LIGASANO® green unsterile	55 x 45 x 2 cm	2
111725176946	16001-007	LIGASANO® green unsterile	55 x 45 x 2 cm	7
111725177512	16001-028	LIGASANO® green unsterile	55 x 45 x 2 cm	28
111725178175	16006-001	LIGASANO® green unsterile	190 x 90 x 2 cm	1
111725179865	16008-006	LIGASANO® green unsterile	190 x 45 x 2 cm	6
111725180659	16110-001	LIGASANO® green unsterile	200 x 11 x 1 cm	1
111725181225	16210-002	LIGASANO® green unsterile	49 x 49 x 9/1 cm	2
111725182915	16812-010	LIGASANO® green sterile	30 x 1.5 x 0.6 cm	10
111725183578	16812-028	LIGASANO® green sterile	30 x 1.5 x 0.6 cm	28
111725184144	16814-010	LIGASANO® green sterile	5 x 5 x 2 cm	10
111725185834	16814-015	LIGASANO® green sterile	5 x 5 x 2 cm	15
111725186400	16814-030	LIGASANO® green sterile	5 x 5 x 2 cm	30
111725187063	16815-010	LIGASANO® green sterile	15 x 5 x 0.6 cm	10
111725188753	16816-010	LIGASANO® green sterile	15 x 10 x 1 cm	10
111725189319	16816-020	LIGASANO® green sterile	15 x 10 x 1 cm	20
111725190113	16816-140	LIGASANO® green sterile	15 x 10 x 1 cm	140
111725191803	16817-010	LIGASANO® green sterile	15 x 10 x 2 cm	10
111725192466	16817-070	LIGASANO® green sterile	15 x 10 x 2 cm	70
111725193032	16820-005	LIGASANO® green sterile	24 x 16 x 2 cm	5
111725194722	16820-035	LIGASANO® green sterile	24 x 16 x 2 cm	35
111725199994	16865-010	LIGASANO® green sterile	200 x 2 x 0.6 cm	10
111725200706	19100-002	LIGASANO® white unsterile	Ø 15 x 8/5 cm	2
111725201369	19100-024	LIGASANO® white unsterile	Ø 15 x 8/5 cm	24
111725203625	36050-006	LIGASANO® Combi pack	5-fold adjustable	1+5
111725204288	36812-010	LIGASANO® orange sterile	30 x 1.5 x 0.6 cm	10
111725205978	36812-028	LIGASANO® orange sterile	30 x 1.5 x 0.6 cm	28
111725207110	36814-010	LIGASANO® orange sterile	5 x 5 x 2 cm	10
111725208800	36814-015	LIGASANO® orange sterile	5 x 5 x 2 cm	15
111725209463	36814-030	LIGASANO® orange sterile	5 x 5 x 2 cm	30
111725212513	36815-005	LIGASANO® Combi pack	5 x 5 x 2 cm	5 x 3
111725213176	36817-010	LIGASANO® orange sterile	15 x 10 x 2 cm	10
111725214866	36820-005	LIGASANO® orange sterile	24 x 16 x 2 cm	5
1117252685871	15072-018	LIGASANO® white unsterile	20 x 10 x 2 cm	18
111767494909	36060-006	LIGASANO® Combi pack	4-fold adjustable	6
111767495572	12013-003	LIGASANO® Combi pack	Various sizes	3
111767496138	12017-022	LIGASANO® Combi pack	Various sizes	22
111804476512	15217-010	LIGASANO® white unsterile	Ø 1.6 x 12 cm	10
111804475946	15218-009	LIGASANO® white unsterile	Ø 2.5 x 12 cm	9

Excerpt from the instructions for use

LIGASANO® white sterile

Please read the complete instructions for use, which is enclosed with each original pack, carefully before use!

General Information:

LIGASANO® white is a therapeutically effective PUR foam material with a wide range of applications for wound care and prevention. The therapeutic uses are based on three underlying physical principles:

Mechanical stimulus: Stimulates local blood circulation within the wound and on the skin, thereby improving the supply of nutrients and oxygen in the wound area (activation of the wound) and acts as a preventative on intact skin. Stimulates normal bodily functions on a local level, especially in the case of inactive patients or, where these functions are reduced. The mechanical stimulus lasts up to three days, by which time LIGASANO® white should then have been changed.

Low-pressure tension: Low-pressure adaption to contours and minimised pressure on wound and body.

Specific suction effect: Absorbs excess fluid without a drying-out effect. Creates a moist and warm wound environment, cleans the wound and thus reduces the bacterial load, without sticking to the wound. Averts skin maceration and is skin-friendly.

Intended Use:

Promotion of wound healing in **sterile application** by **cleaning, covering or filling wounds** or in non-sterile application as a secondary dressing (without direct wound contact) in wound treatment or in concomitant wound treatment, e.g. for pressure relief, or preventively as pressure and friction protection.

Indications:

Acute and chronic wounds (e.g. decubital ulcers, ulcer cruris, diabetic foot ulcers). Post-surgical wounds and post-surgical wound healing impairments. Thermal wounds, such as burns, scalds, frostbite

Adverse Effects, Side Effects, Residual Risks:

Wound cleaning: mechanical debridement may cause pain, even with LIGASANO® white sterile.

Wound treatment: because of its cell structure a vascularisation (e.g. of granulation tissue) into the material is possible. You can prevent this with a timely dressing change, after 1-3 days using LIGASANO® white sterile. On passive wounds or parts of the body that receive poor blood supply, pain may occur after application. This occurs when the blood circulation is stimulated by the mechanical stimulus of the material to such an extent, that previously diminished or reduced pain returns. Initially the pain can be intense, however it is usually regulated after a few hours to days. If this effect is not desired, you should forgo its application. Too much pressure due to the material's prior tension, external influence, or suction may lead to compression of vessels and therefore to pressure ulcers.

It may occur that LIGASANO® white sterile sticks to the wound, mostly due to too much air reaching the wound. This means that the LIGASANO® white sterile dressing is not thick enough (it should be 2 cm from skin level). LIGASANO® white can be moistened on the side facing the wound (e.g. with Ringer's solution or wound irrigation solution) to prevent it from sticking.

If hypergranulation occurs, please treat it appropriately and switch to other dressing material groups, if necessary.

LIGASANO® white is a very simple product and its effects are due to its structure. It does not contain nor emit any active compounds. With the correct application of LIGASANO® white there are no known negative side effects, nor incompatibility or interactions with medicine. Typical side effects may include reddening of the skin upon contact and initial "itching". More intense side effects may include return of sensitivity due to stimulation of local blood circulation. Its effects are only physical and end almost immediately after discontinuation of application. As with all wound dressings, minor skin reactions such as maceration, erythema, secondary infected dermatitis, erysipelas as well as hypersensitivity reactions and pain at dressing change may occur.

Please report any incidents that have occurred in connection with our dressing material LIGASANO® white sterile to us or to the competent authority of the member state in which the incident occurred.

Contraindications:

Tumour wounds and untreated osteomyelitis. Not suitable for contact with organs, contact to exposed blood vessels with the danger of damage, or contact with nerves.

Patient Target Group:

Children (> 2 - 16 years of age) and adults (> 16 years of age) with acute, chronic, post-surgical wounds/wound healing disturbances and thermal wounds.

User/Practitioner:

Application must always be done by healthcare professionals. The decision about the particular kind of application of LIGASANO® white sterile is determined by healthcare professionals and is adapted and appropriated in every single case.

Capability Characteristics and Clinical Benefit:

With the sterile wound dressing LIGASANO® white you can clean, cover or pack wounds. The product can be used both for mechanical debridement and for wound treatment in all phases of wound healing. The expected clinical benefit presents itself in a very effective wound cleaning, wound bed conditioning and in the healing of the wounds.

LIGASANO® white unsterile

Please read the complete instructions for use, which is enclosed with each original pack, carefully before use!

General Information:

LIGASANO® white is a therapeutically effective PUR foam material with a wide range of applications for wound care and prevention. The therapeutic uses are based on three underlying physical principles:

Mechanical stimulus: Stimulates local blood circulation within the wound and on the skin, thereby improving the supply of nutrients and oxygen in the wound area (activation of the wound) and acts as a preventative on intact skin. Stimulates normal bodily functions on a local level, especially in the case of inactive patients or, where these functions are reduced. The mechanical stimulus lasts up to three days, by which time LIGASANO® white should then have been changed.

Low-pressure tension: Low-pressure adaption to contours and minimised pressure on wound and body.

Specific suction effect: Absorbs excess fluid without a drying-out effect. Creates a moist and warm wound environment, cleans the wound and thus reduces the bacterial load, without sticking to the wound. Averts skin maceration and is skin-friendly.

Intended Use:

In **non-sterile application** as a **secondary dressing (without direct wound contact)** in wound treatment or in **concomitant wound treatment**, e.g. for pressure relief, or **preventively** as **pressure and friction protection**.

Indications:

Concomitant wound therapy in case of acute and chronic wounds (e.g. decubital ulcers, ulcer cruris, diabetic foot ulcers). Also post-surgical wounds and post-surgical wound healing impairments and thermal wounds, such as burns, scalds, frostbite. Preventive use when the skin is still intact, e.g. pressure ulcer prophylaxis through active positioning on LIGASANO® white unsterile. Prophylaxis of maceration and intertrigo.

Adverse Effects, Side Effects, Residual Risks

There are no particular risks when using LIGASANO® white unsterile for pressure relief and hygiene, if the general nursing rules and our instructions are noted. LIGASANO® white unsterile is a very simple product and its effects are due to its structure. It does not contain nor emit any active compounds. With the correct application of LIGASANO® white there are no known negative side effects, nor incompatibility or interactions with medicine. Typical side effects may include reddening of the skin upon contact and initial "itching". Its effects are only physical and end almost immediately after discontinuation of application.

Please report any incidents that have occurred in connection with our dressing material LIGASANO® white unsterile to us or to the competent authority of the member state in which the incident occurred.

Contraindications:

Very sensible skin (e.g. parchment skin), direct contact to wounds.

Patient Target Group:

Children (> 2 - 16 years of age) and adults (> 16 years of age) with or without wounds.

User/Practitioner:

Application must always be done by healthcare professionals. The decision about the particular kind of application of LIGASANO® white unsterile is determined by healthcare professionals and is adapted and appropriated in every single case.

Capability Characteristics and Clinical Benefit:

With the unsterile wound dressing LIGASANO® white you can prevent the risk of the development of pressure ulcer and moisture-caused wounds. It is able to prevent the development of intertrigo and partly fungal infections, too. The expected clinical benefit present itself in a very effective pressure relief. Both in already existing wounds and preventive. An additional benefit is the moisture-regulating effect.

Excerpt from the instructions for use

LIGASANO® green sterile

Please read the complete instructions for use, which is enclosed with each original pack, carefully before use!

General Information:

LIGASANO® green is an elastic, expanded polyurethane (PUR), which is almost always pervious to air and liquids (e.g. water) - imagine a three-dimensional sieve. It has a coarse texture and rough surface.

Intended Use:

Promotion of wound healing in **sterile application** by **cleaning, covering or filling wounds** or in non-sterile application as a secondary dressing (without direct wound contact) in wound treatment or in concomitant wound treatment, e.g. for pressure relief, or preventively as pressure and friction protection.

Indications:

Acute and chronic wounds (e.g. decubital ulcers, ulcus cruris, diabetic foot ulcers). Post-surgical wounds and post-surgical wound healing impairments. Thermal wounds, such as burns, scalds, frostbite.

Adverse Effects, Side Effects, Residual Risks:

Wound cleaning: mechanical debridement may cause pain, even with LIGASANO® green sterile.

Wound treatment: because of its cell structure a vascularisation (e.g. of granulation tissue) into the material is possible. You can prevent this with a timely dressing change, after 1-4 days using LIGASANO® green sterile. On passive wounds or parts of the body that receive poor blood supply, pain may occur after application. This occurs when the blood circulation is stimulated by the mechanical stimulus of the material to such an extent, that previously diminished or reduced pain returns. Initially the pain can be intense, however it is usually regulated after a few hours to days. If this effect is not desired, you should forgo its application. In case of an "open wound care" without additional cover you have to take into consideration that there exists barely protection against the ambient atmosphere. Too much pressure due to the material's prior tension, external influence, or suction may lead to compression of vessels and therefore to pressure ulcers.

If hypergranulation occurs, please treat it appropriately and switch to other dressing material groups, if necessary.

LIGASANO® green is a very simple product and its effects are due to its structure. It does not contain nor emit any active compounds. With the correct application of LIGASANO® green there are no known negative side effects, nor incompatibility or interactions with medicine. Typical side effects may include reddening of the skin upon contact and initial "itching". More intense side effects may include return of sensitivity due to stimulation of local blood circulation. Its effects are only physical and end almost immediately after discontinuation of application. As with all wound dressings, minor skin reactions such as maceration, erythema, secondary infected dermatitis, erysipelas as well as hypersensitivity reactions and pain at dressing change may occur.

Please report any incidents that have occurred in connection with our dressing material LIGASANO® green sterile to us or to the competent authority of the member state in which the incident occurred.

Contraindications:

Tumour wounds and untreated osteomyelitis. Not suitable for contact with organs, contact to exposed blood vessels or anastomotic areas with the danger of damage, or contact with nerves.

Patient Target Group:

Children (> 2 -16 years of age) and adults (> 16 years of age) with acute, chronic, post-surgical wounds/wound healing disturbances and thermal wounds.

User/Practitioner:

Application must always be done by healthcare professionals. The decision about the particular kind of application of LIGASANO® green sterile is determined by healthcare professionals and is adapted and appropriated in every single case.

Capability Characteristics and Clinical Benefit:

With the sterile wound dressing LIGASANO® green you can **clean, cover or pack wounds**. The product can be used both for **mechanical debridement in all phases of wound healing** and for **wound treatment during the phase of exudation and granulation**. The expected clinical benefit present itself in a very effective wound cleaning and wound bed conditioning.

LIGASANO® green unsterile

Please read the complete instructions for use, which is enclosed with each original pack, carefully before use!

General Information:

LIGASANO® green is an elastic, expanded polyurethane (PUR), which is almost always pervious to air and liquids (e.g. water) - imagine a three-dimensional sieve. It has a coarse texture and rough surface.

Intended Use:

In **non-sterile application** as a **secondary dressing (without direct wound contact)** in wound treatment or in **concomitant wound treatment**, e.g. for pressure relief, or **preventively as pressure and friction protection**.

Indications:

Concomitant wound therapy in case of acute and chronic wounds (e.g. decubital ulcers, ulcus cruris, diabetic foot ulcers). Also post-surgical wounds and post-surgical wound healing impairments and thermal wounds, such as burns, scalds, frostbite. Preventive use when the skin is still intact, e.g. pressure ulcer prophylaxis through active positioning on LIGASANO® green unsterile. Prophylaxis of maceration.

Adverse Effects, Side Effects, Residual Risks:

There are no particular risks when using LIGASANO® green unsterile for pressure relief and hygiene, if the general nursing rules and our instructions are noted.

LIGASANO® green unsterile is a very simple product and its effects are due to its structure. It does not contain nor emit any active compounds. With the correct application of LIGASANO® green there are no known negative side effects, nor incompatibility or interactions with medicine.

Please report any incidents that have occurred in connection with our dressing material LIGASANO® green unsterile to us or to the competent authority of the member state in which the incident occurred.

Contraindications:

Direct contact to skin, direct contact to wounds

Patient Target Group:

Children (> 2 - 16 years of age) and adults (> 16 years of age) with or without wounds.

User/Practitioner:

Application must always be done by healthcare professionals. The decision about the particular kind of application of LIGASANO® green unsterile is determined by healthcare professionals and is adapted and appropriated in every single case.

Capability Characteristics and Clinical Benefit:

With the unsterile wound dressing LIGASANO® green you can prevent the risk of the development of pressure ulcer and moisture-caused wounds.

The expected clinical benefit present itself in a very effective pressure relief. Both in already existing wounds and preventive. An additional benefit is the moisture-regulating effect.

Excerpt from the instructions for use

LIGASANO® orange sterile

Please read the complete instructions for use, which is enclosed with each original pack, carefully before use!

General Information:

LIGASANO® orange is an elastic, expanded polyurethane (PUR), which is almost always pervious to air and liquids (e.g. water) - imagine a three-dimensional sieve. It has a coarse texture and a relatively rough surface.

Intended Use:

Promotion of wound healing in **sterile application** by **cleaning, covering or filling wounds** or in non-sterile application as a secondary dressing (without direct wound contact) in wound treatment or in concomitant wound treatment, e.g. for pressure relief, or preventively as pressure and friction protection.

Indications:

Acute and chronic wounds (e.g. decubital ulcers, ulcus cruris, diabetic foot ulcers). Post-surgical wounds and post-surgical wound healing impairments. Thermal wounds, such as burns, scalds, frostbite.

Adverse Effects, Side Effects, Residual Risks:

Wound cleaning: mechanical debridement may cause pain, even with LIGASANO® orange sterile.

Wound treatment: because of its cell structure a vascularisation (e.g. of granulation tissue) into the material is possible. You can prevent this with a timely dressing change, after 1-4 days using LIGASANO® orange sterile. On passive wounds or parts of the body that receive poor blood supply, pain may occur after application. This occurs when the blood circulation is stimulated by the mechanical stimulus of the material to such an extent, that previously diminished or reduced pain returns. Initially the pain can be intense, however it is usually regulated after a few hours to days. If this effect is not desired, you should forgo its application. In case of an "open wound care" without additional cover you have to take into consideration that there exists barely protection against the ambient atmosphere. Too much pressure due to the material's prior tension, external influence, or suction may lead to compression of vessels and therefore to pressure ulcers.

If hypergranulation occurs, please treat it appropriately and switch to other dressing material groups, if necessary.

LIGASANO® orange is a very simple product and its effects are due to its structure. It does not contain nor emit any active compounds. With the correct application of LIGASANO® orange there are no known negative side effects, nor incompatibility or interactions with medicine.

Typical side effects may include reddening of the skin upon contact and initial "itching". More intense side effects may include return of sensitivity due to stimulation of local blood circulation. Its effects are only physical and end almost immediately after discontinuation of application.

Please report any incidents that have occurred in connection with our dressing material LIGASANO® orange sterile to us or to the competent authority of the member state in which the incident occurred.

Contraindications:

Tumour wounds and untreated osteomyelitis. Not suitable for contact with organs, contact to exposed blood vessels or anastomotic areas with the danger of damage, or contact with nerves.

Patient Target Group:

Children (> 2 - 16 years of age) and adults (> 16 years of age) with acute, chronic, post-surgical wounds/wound healing disturbances and thermal wounds.

User/Practitioner:

Application must always be done by healthcare professionals. The decision about the particular kind of application of LIGASANO® orange sterile is determined by healthcare professionals and is adapted and appropriated in every single case.

Capability Characteristics and Clinical Benefit:

With the sterile wound dressing LIGASANO® orange you can **clean, cover or pack wounds**. The product can be used both for **mechanical debridement in all phases of wound healing** and for **wound treatment during the phase of exudation and granulation**. The expected clinical benefit present itself in a very effective wound cleaning and wound bed conditioning.

LIGASANO® orange unsterile

Please read the complete instructions for use, which is enclosed with each original pack, carefully before use!

General Information:

LIGASANO® orange is an elastic, expanded polyurethane (PUR), which is almost completely pervious to air and liquids - imagine a three-dimensional sieve. It has a coarse texture and a relatively rough surface. LIGASANO® orange unsterile is used for concomitant wound treatment and care, if a strength and elastic material with high permeability to air and liquid is required

Intended Use:

In **non-sterile application** as a **secondary dressing (without direct wound contact)** in wound treatment or in **concomitant wound treatment**, e.g. for pressure relief, or **preventively as pressure and friction protection**.

Indications:

Concomitant wound therapy in case of acute and chronic wounds (e.g. decubital ulcers, ulcus cruris, diabetic foot ulcers). Also post-surgical wounds and post-surgical wound healing impairments and thermal wounds, such as burns, scalds, frostbite.

Preventive use when the skin is still intact, e.g. pressure ulcer prophylaxis through active positioning on LIGASANO® orange unsterile. Prophylaxis of maceration.

Adverse Effects, Side Effects, Residual Risks:

There are no particular risks when using LIGASANO® orange unsterile for pressure relief and hygiene, if the general nursing rules and our instructions are noted.

LIGASANO® orange unsterile is a very simple product and its effects are due to its structure. It does not contain nor emit any active compounds. With the correct application of LIGASANO® orange unsterile there are no known negative side effects, nor incompatibility or interactions with medicine.

Please report any incidents that have occurred in connection with our dressing material LIGASANO® orange unsterile to us or to the competent authority of the member state in which the incident occurred.

Contraindications:

Sensible skin (e.g. parchment skin), direct contact to wounds.

Patient Target Group:

Children (> 2 - 16 years of age) and adults (> 16 years of age) with or without wounds.

User/Practitioner:

Application must always be done by healthcare professionals. The decision about the particular kind of application of LIGASANO® orange unsterile is determined by healthcare professionals and is adapted and appropriated in every single case.

Capability Characteristics and Clinical Benefit:

With the unsterile wound dressing LIGASANO® orange you can prevent the risk of the development of pressure ulcer and moisture-caused wounds.

The expected clinical benefit present itself in a very effective pressure relief. Both in already existing wounds and preventive. An additional benefit is the moisture-regulating effect.

Prices are subject to change

Terms of delivery:

Ex work Cadolzburg/Germany or on appointment.

Terms of payment:

Payment in advance or on appointment.

Place of jurisdiction:

Fuerth or Nuremberg / Germany